

# Tranquility® Premium TopLiner® Booster Pads

Flow-through pad designed to add capacity to any disposable incontinence product. Thin design allows for placement in target areas where additional absorption is needed. Secures with adhesive strip.



premium

| Item No. | Size    | Length/Width  | Capacity       | Case Count         | HCPC Code | Best Suited for Host Garment Sizes |
|----------|---------|---------------|----------------|--------------------|-----------|------------------------------------|
| 2072     | Mini    | 10.5" x 2.75" | 5.7 oz/169 mL  | 200 (8 bags of 25) | T4535     | Youth, XS                          |
| 2070     | Regular | 14" x 4"      | 10.9 oz/322 mL | 200 (8 bags of 25) | T4535     | S, M, L                            |
| 2060     | Super   | 15" x 4.25"   | 14.2 oz/420 mL | 200 (8 bags of 25) | T4535     | XL, 2XL                            |



ABSORBENCY:  
14.2 oz/420 mL (Size Regular)



Booster pads have a flow-through design (no moisture barrier) that allows the pad to absorb initial fluid, then pass additional fluid through to the primary undergarment. They must be used in combination with another disposable incontinence product that has a moisture-proof backing and should not be used with regular underwear.

### Benefits of Booster Pads

- Adds capacity – Extends weartime allowing for uninterrupted sleep and social activities
- Discreet – Clothlike material is gentle against the skin and quiet during movement.
- Helps prevent leakage from top, sides and leg openings
- Reduces changes of the disposable undergarment
- Secure placement
- Not made with natural rubber latex

### Features

Flow-Through Design • Clothlike Adhesive Strip • Embossed Channels





Want more information? 1-800-467-3224  
TranquilityProducts.com



# What is a Booster Pad?

**Booster pads have a flow-through design** allowing them to absorb initial fluid before passing additional moisture into a primary incontinence product. Booster pads should always be placed inside another disposable incontinence product that has a moisture-proof backing. **Do not use with regular underwear.**

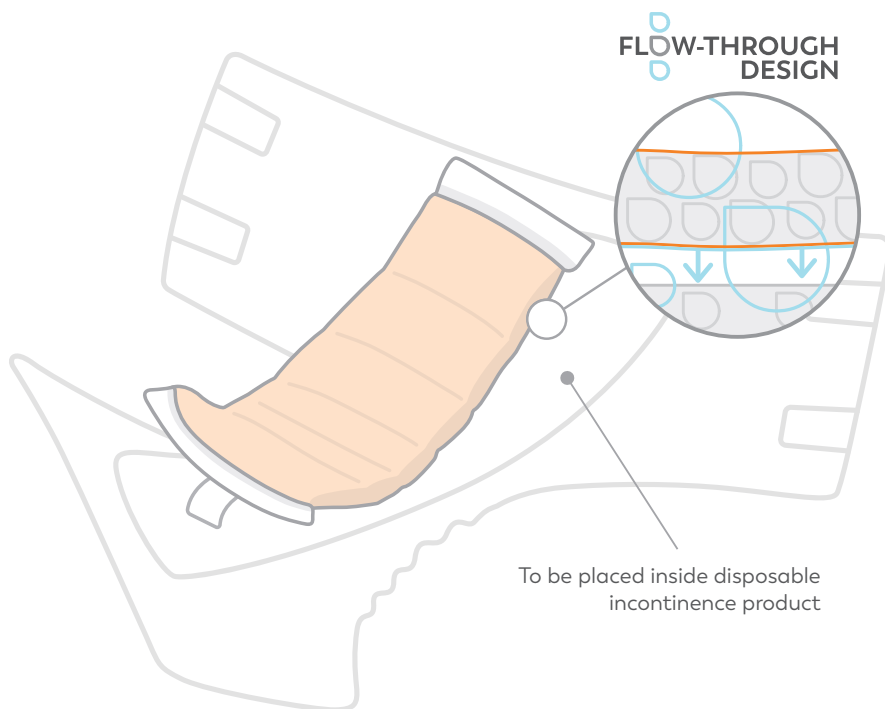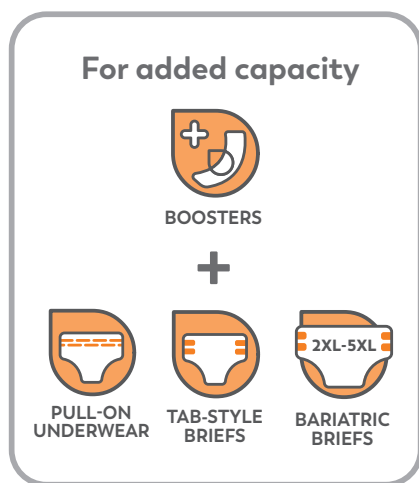

## Did You Know?

Booster pad placement is important. **The pads can be placed in targeted areas where additional absorption is needed.**



Booster pads can be folded and even doubled for added protection, but should never be cut.



**Host Garment** – The pad's adhesive side should face out, away from the body. Make sure the pad is fully contained within the host incontinence product; no part should extend beyond the edges.

**For Women** – Place the pad more toward the bottom.

**For Men** – Place the pad more toward the front.



Want more information?  
1-800-467-3224 • [TranquilityProducts.com](https://TranquilityProducts.com)

